


# Easy Catch<sup>+</sup>

HAEMOSTATIC Y-ADAPTOR




Easy, Fast, Safe & Accurate

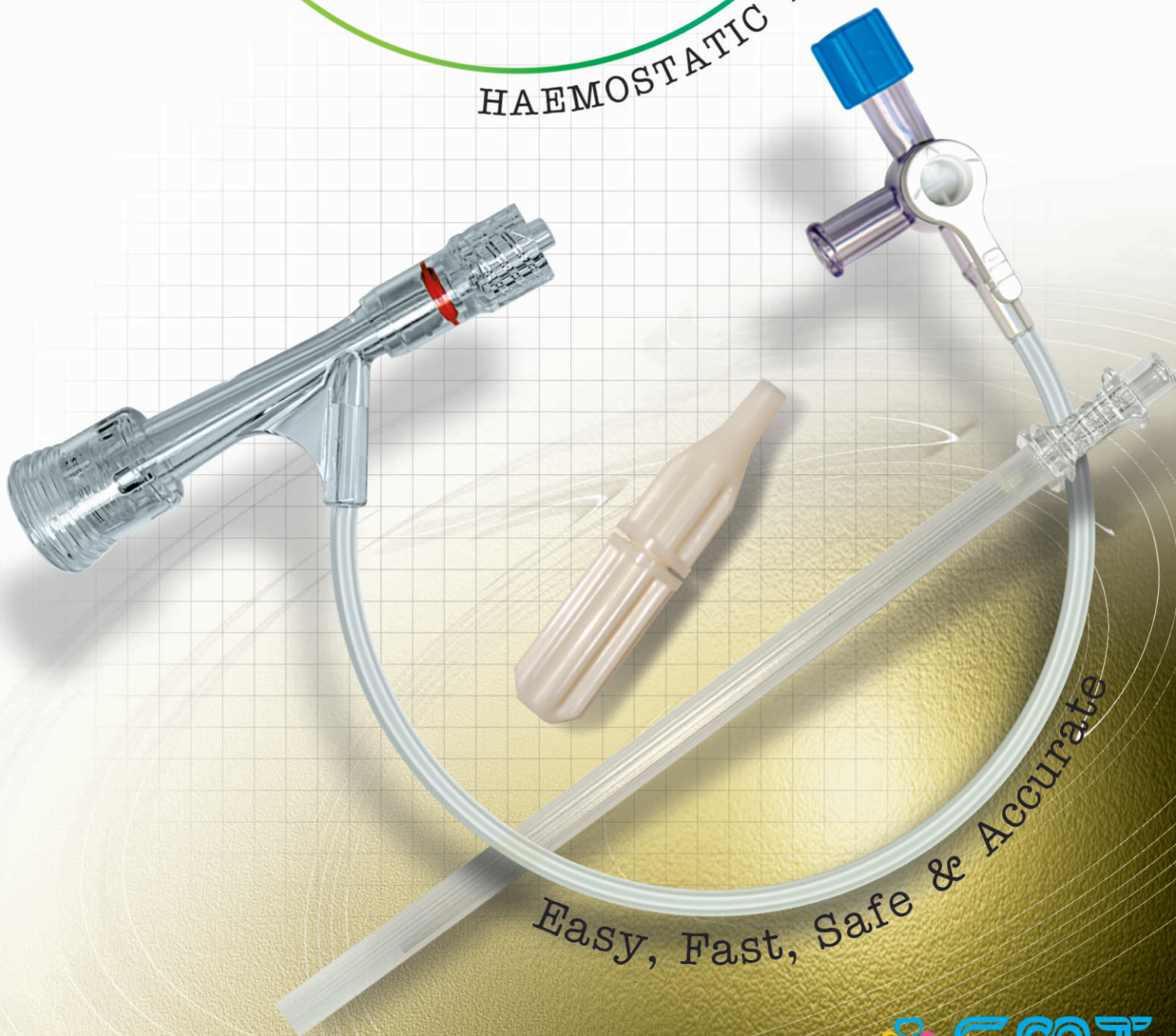

# Easy Catch<sup>+</sup>

## HAEMOSTATIC Y-ADAPTOR

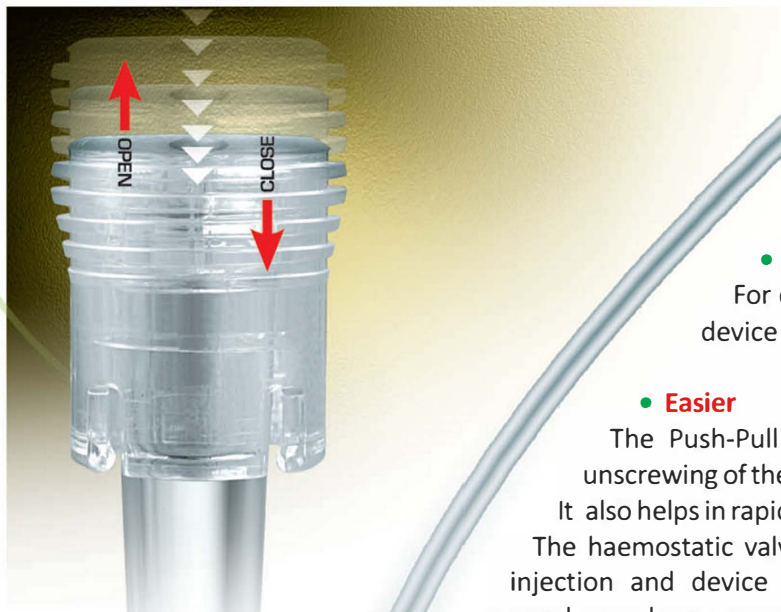

The Easy-Catch+ with its unique Push-Pull button valve opener is easier, safer & cleaner to use.

- **A Unique design**

For easy, fast and accurate single handed intra vascular device insertion.

- **Easier**

The Push-Pull mechanism eliminates the need of screwing & unscrewing of the valve.

It also helps in rapid and efficient device insertion.

The haemostatic valve system permits simultaneous contrast medium injection and device positioning, a major advantage during stenting procedures where accurate placement is required.

- **Essential for DES**

Perfect insertion device respecting the integrity of the coated surface.

- **Clean and Safer**

Combination of a haemostatic valve and an integrated Push-Pull button valve opener. Cleaner and safer procedure with a minimized backflow thus limiting blood loss.

Y-Connector	
Inner Lumen Compatibility	7F / 0.092" / 2.33 mm maximum
Pressure Resistance (1) after introduction of a device	8 atm. maximum
Metallic Insertion Tool	
External Diameter	2.7 F / 0.90 mm

Inner Lumen Diameter	1.9 F / 0.025" / 0.63 mm
Usable Length	10 cm
Torquer	
Guidewire Compatibility	From 0.010" minimum to 0.035" maximum

*\*Resistance of the haemostatic valve to the pressure of injection through the sideport*



Sahajanand Medical Technologies Pvt. Ltd.  
 "Sahajanand Estate", Wakharia Wadi,  
 Near Dabholi Char Rasta,  
 Ved Road, Surat 395004, Gujarat, INDIA  
 Tel.: +91 261 3060606 Fax: +91 261 3060607  
 Email : contact@sahmed.com  
 Visit us at : www.smtpl.com  
 • CIN: U33119GJ2001PTC040121