




MeCross CTO ϕ 2.0-4.0mm illustration

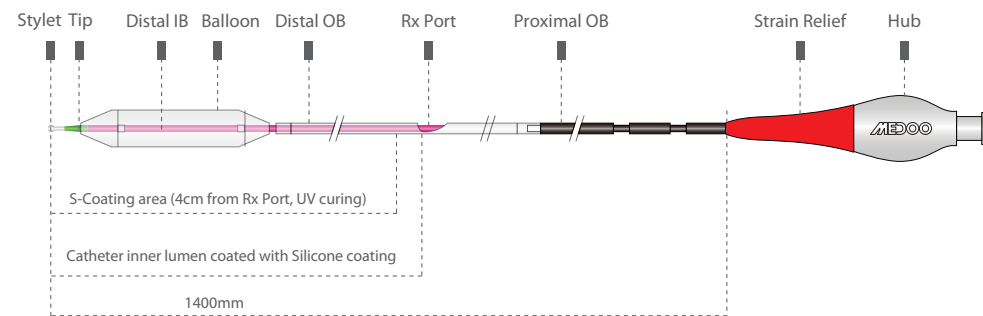

MeCross CTO ϕ 1.0-1.75mm distal illustration

ORDERING INFORMATION

ϕ (mm)	1.0	1.25	1.5	1.75	2.0	2.25	2.5	2.75	3.0	3.25	3.5	4.0
Length(mm)												
10	MC-01-1010	-1210	-1510	-1710	-2010	-2210	-2510	-2710	-3010	-3210	-3510	-4010
15	MC-01-1015	-1215	-1515	-1715	-2015	-2215	-2515	-2715	-3015	-3215	-3515	-4015
20	MC-01-1020	-1220	-1520	-1720	-2020	-2220	-2520	-2720	-3020	-3220	-3520	-4020
30	MC-01-1030	-1230	-1530	-1730	-2030	-2230	-2530	-2730	-3030	-3230	-3530	-4030

COMPLIANCE INFORMATION

ϕ (mm)	1.0	1.25	1.5	1.75	2.0	2.25	2.5	2.75	3.0	3.25	3.5	4.0
Pressure(atm)												
2	0.90	1.15	1.40	1.65	1.87	2.15	2.36	2.60	2.82	3.05	3.30	3.78
4	0.95	1.20	1.45	1.70	1.94	2.20	2.45	2.67	2.91	3.15	3.40	3.90
6(NP)	1.00	1.25	1.50	1.75	2.00	2.25	2.50	2.75	3.00	3.25	3.50	4.00
8	1.05	1.30	1.55	1.80	2.05	2.35	2.57	2.85	3.10	3.35	3.63	4.13
10	1.10	1.35	1.60	1.85	2.11	2.40	2.64	2.95	3.18	3.45	3.74	4.25
12	1.15	1.40	1.65	1.90	2.19	2.50	2.71	3.00	3.27	3.55	3.85	4.35
14	1.20	1.45	1.70	1.95	2.25	2.55	2.80	3.05	3.35	3.65	3.95	4.45
16(RBP)	1.23	1.50	1.75	2.00	2.32	2.60	2.85	3.15	3.45	3.75	4.05	4.58
18	1.25	1.55	1.80	2.05	2.38	2.68	2.93	3.25	3.55	3.85	4.15	4.68



About Medoo

Medoo Medical is a high-tech enterprise which focuses on R&D, production and sales of advanced intravascular/extravascular intervention-al medical devices. Founded in 2014, the Medoo group consists of Shenzhen Medoo Medical Tech. Co., Ltd., Shenzhen Medoo Medtech. Co., Ltd., Jiangsu Medoo Medical Equipment Tech. Co., Ltd. and Hong Kong Medoo Medical Tech. Co., Ltd.

Medoo's team boasts solid technical capabilities and strong policy support from the Shenzhen Science and Technology Commission. Centered on cardiovascular and peripheral vascular intervention fields, Medoo has developed many core products: MeCross CTO /NC /HP dilatation balloon catheter, MeCross PTA peripheral balloon catheter, NeuCross RX /OTW intracranial balloon catheter, MePower PTA peripheral high-pressure balloon catheter, arterial compression tourniquet, guidewire, zebra guidewire, bronchial occlusion catheter, Lebon combo coronary balloon, SkyRH leading micro-catheter, single-wire/dual-wire scoring balloon.

Integrating the resources of clinical experts in cardiovascular intervention, scientists of Chinese Academy of Science and the R&D specialists on medical consumables of intravascular intervention, Medoo Medical has launched industry-university-research cooperation united with numerous universities and hospitals and established international-standard R&D and production centers in Shenzhen and Jiangsu respectively. Medoo Medical has obtained BSI ISO 13485 certification (MD: 662153) and GMP certification of China. Its flagship products deliver reliable performance and hold EU CE certification and Chinese NMPA certification, earning recognition from global clients. Medoo Medical not only carries out distribution business in the world with its own brand, but also provides professional OEM/ODM/CDMO production technical services for high-end vascular intervention medical device enterprises at home and abroad.

Culture of Medoo

Craftsmanship Spirit

Pursue excellence in technology, strive for perfection and develop top level products.

Commitment to Quality

Strictly abide to laws and regulations, in order to ensure the health & safety of patients

Innovation for the Future

Focus on the future life science and technology, create a vigorous Medoo brand

Treasure the Life

Improve patients' physical & mental health, care for human beings, contribute to our society

Goals of Medoo

- Develop Medoo Medical to the world's leading national high-end medical device brand
- Benefiting millions of patients with vascular diseases
- Pursue the harmony and unity of the enterprise, patients and society

Shenzhen Medoo Medical Tech. Co., Ltd.

Room A-408 , Platform Building, Shenzhen Virtual University Park, Yuexing 2 Road, Nanshan district, Shenzhen, China

+86 755-26928966

www.medoo.hk

Jiangsu Medoo Medical Equipment Tech. Co., Ltd.

NO.32 Century Avenue, Dongkan Town, Binhai County, Yancheng, Jiangsu, 224500, China

+86 515-84537776

zxy@medoo.hk

Semi-Compliant PTCA Balloon Dilatation Catheter

NMPA: 国械注准20183030557

CE662148

MDR 786405

"Seamless" Force Conduction
supreme pushability

Superior Hydrophilic "S-Coating"

ensures durable and strong hydrophilic ability to pass through the most complex lesions

Minimum Diameter of 1.0mm

designed for chronic total occlusion stenosis

Delicate "A-Advance" Tip

excellent crossing performance






Semi-Compliant PTCA Balloon Dilatation Catheter [MeCross[®] CTO]

MeCross CTO (RX type) Semi-Compliant PTCA Balloon Dilatation Catheter (NMPA: 国械注准 20183030557, CE662148, MDR786405) is a core product of Medoo Medical MeCross[®] series, it is designed for CTO and other chronic vascular stenotic lesions. The latest technologies in balloon design and materials are applied to ensure excellent crossing performance and stable intraoperative safety.

MeCross CTO applies a thin OD pusher shaft design, enabling two balloons to be delivered simultaneously inside one 2 mm (6F) guiding catheter for kissing dilatation surgery.

The unique A-Advance Tip design ensures the product with superior crossing ability, while the new ultra-lubricated hydrophilic S-Coating also helps reduce shaft friction and further enhance lesion crossing performance in PTCA surgery. "Seamless" force conduction technology ensures the whole catheter a satisfied balance in the pushing and tracking ability, which makes MeCross CTO able to pass through stenosis smoothly and deal with the hard calcification lesion in CTO case.

The diameter specifications which ranges from 1.0mm to 4.0mm meets most pre-dilatation clinical needs of cardiovascular interventional surgery.

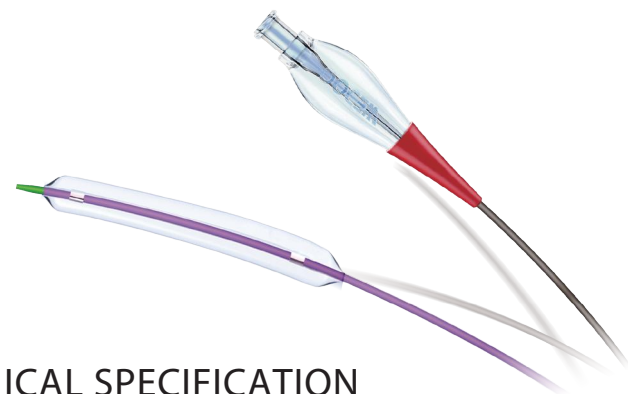

TECHNICAL SPECIFICATION

RX Design

Working length 140cm


12 diameter configurations 1.0mm-4.0mm

4 length configurations 10mm, 15mm, 20mm, 30mm


NP=6atm, RBP=16atm

Product Advantages


- Specially designed for CTO (1.0-1.25mm) and other chronic coronary stenosis pre-dilatation case
- Advantage of optimized CTO case configurations (1.0-1.25mm)
 - Suitable to extremely challenging lesions with minimum diameter of 1.0mm
 - Single marker design and laser rotary cut distal core attribute to excellent trackability and flexibility
 - Handle the most complex CTO lesion, unlimited crossing force
- "Seamless" force conduction technology ensures enough pushability
- "S-Coating" with the latest technology can ensure durable and strong hydrophilic ability to pass through the most complex lesions with zero resistance
- Delicate "A-advance" tip design to make the catheter softer and more flexible to protect the blood vessels and defend the health
- Outstanding kissing capability to 6F guiding catheter

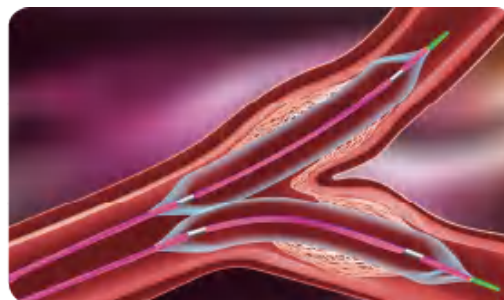


Outstanding kissing capability is suitable for 6F guiding catheter.

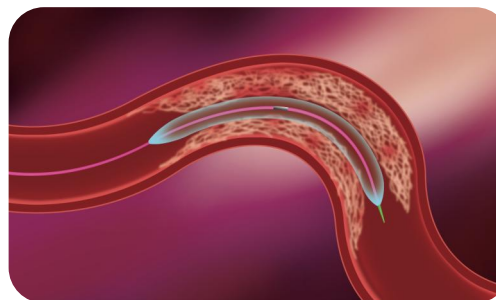

It can replace microcatheter to treat CTO case.

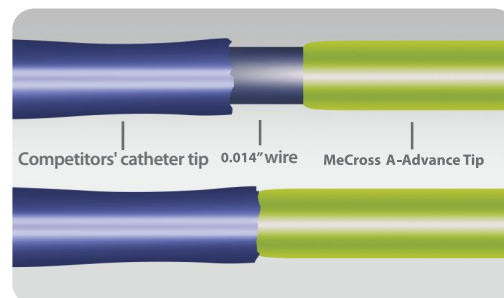

A-Advance Tip successfully traverses the competitors' tip.

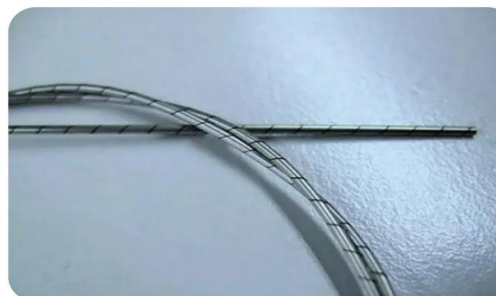

The laser rotary cut distal core keeps good flexibility.


Description of S-Coating

Performance & Advantages

The S-Coating hydrophilic coating adopts the latest UV light curing technology to significantly reduce surface friction and help decrease the trauma to patient. While effectively reducing the patient's discomfort the complications caused by mucosal and tissue cell damages during product insertion are also largely avoided. Moreover, S-Coating can also speed up the delivery of the catheter.

Medoo S-Coating vs. coating of competitor's catheter

Coefficient of Friction(COF)

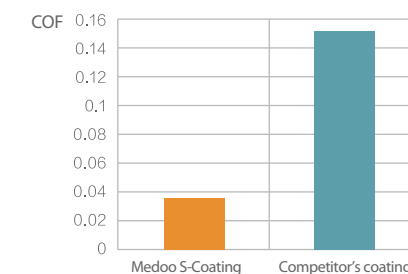

Comparison with Other Competitors

Item	Coating Performance
COF	COF of Medoo S-Coating=0.0353 COF of competitor's coating=0.151
Firmness	After 25 friction test cycles (gripping force=300g), Medoo S-Coating keeps stable and integrated, while competitor's coating is abraded.

Instruction

